

Installation PurifiX® 8

Avant toute utilisation, veuillez lire attentivement les documents liés à l'appareil.

Téléchargez les documents via le **code QR**, ou sur : www.loopdeescience.com/downloads



Réception

- Retirer les films de protection et les adhésifs de maintien de l'appareil.


Bride de transport

Avant toute mise sous tension

- Sous l'appareil, retirer la vis de son **emplacement de transport**, et la placer dans l'**emplacement de stockage** (voir schémas).

Pour transporter l'appareil

- Aligner l'axe des tiges aimantées et consommables peignes à la main, et replacer la vis dans son **emplacement de transport** (voir schémas).

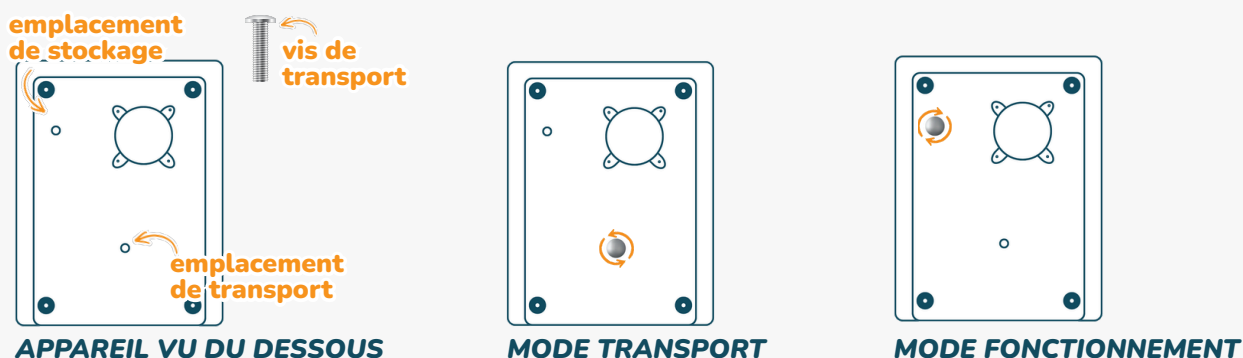

Nettoyage

- A l'aide d'un tissu absorbant imbibé d'alcool ménager (70°), nettoyer les surfaces extérieures, intérieures et le support du PurifiX 8.

Mise sous tension

- Mettre l'appareil sous tension en respectant les indications électriques données dans la notice.

Stérilisation UV

- A l'aide de l'interface de pilotage, menu « **Settings** » et « **UV stérilisation** » : Lancer un cycle de stérilisation de **30 minutes**.

Le PurifiX 8 est prêt à être utilisé selon le protocole défini

Quick start PurifiX® 8

Before using the device, please read the accompanying documentation carefully.

Download the documents via the QR code, or at:
www.loopdeescience.com/downloads



Unpacking:

- Remove the protective films and adhesive strips from the device.


Transport bracket

Before switching on the device:

- Underneath the device, remove the screw from its **transport position** and place it in the **storage position** (see *diagrams*).

To transport the device

- Align the axes of the magnetic rods and comb consumables by hand and replace the screw in its **transport position** (see *diagrams*).

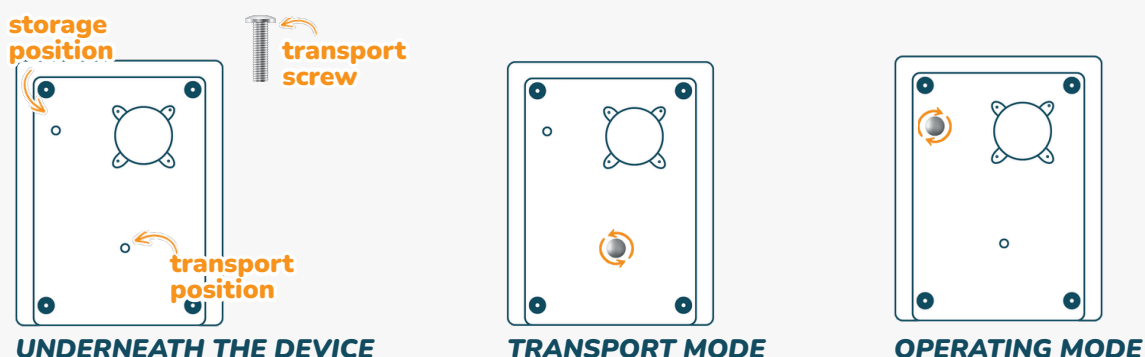

Cleaning

- Using an absorbent cloth soaked in household alcohol (70%), clean the exterior and interior surfaces and the base of the PurifiX 8.

Switching on

- Switch on the device in accordance with the electrical instructions provided in the user manual.

UV Sterilisation

- Using the control interface, select the 'Settings' menu and 'UV Sterilisation': Start a 30-minute sterilisation cycle.

PurifiX 8 is ready for use in accordance with the specified protocol